

TITANIUM Adjustable Total


CODE
TAP 07.63

TI: Titanium ASTM F67 Gr. 2

- WINDOWED (OPENED) FLANGE – it allows the clear view of the below surgery field
- TITANIUM STRUCTURE – Titanium is a bioinert material, largely tested in surgery fields and with high mechanical resistance and lightness
- EXTENSIBLE (ADJUSTABLE) STEM – the stem is extensible through the traction and shortable through the compression
- GOOD SOUND WAVE CONDUCTIVITY – oscillometric tests using Laser/Doppler show the unchanged conductivity of the sound wave in all stem configurations



MRI compatible till 3 Tesla.

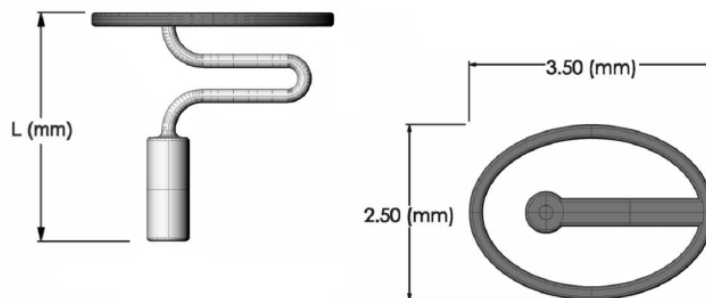

Code	L (mm)
TAP 07.63	From 3.00 to 8.00

Wire \varnothing : 0.20 mm
Columella \varnothing : 0.60 mm
Weight: 5.20 mg

Packaging:
1 product each box