

Audio-HA® Offset Total Prosthesis

CODE
TPL 07.05

HA: Bioactive calcium hydroxylapatite ceramic


It replaces the entire ossicular chain, from the malleus or the tympanic membrane to the stapes footplate.

The interposition of cartilage or bone paste between the flange and the tympanic membrane is not required.

The offset columella remains perpendicular to the oval window.

The shape of the flange makes it possible to see the distal end of the columella.

The flange can be bound to the tympanic membrane by fibrin glue.

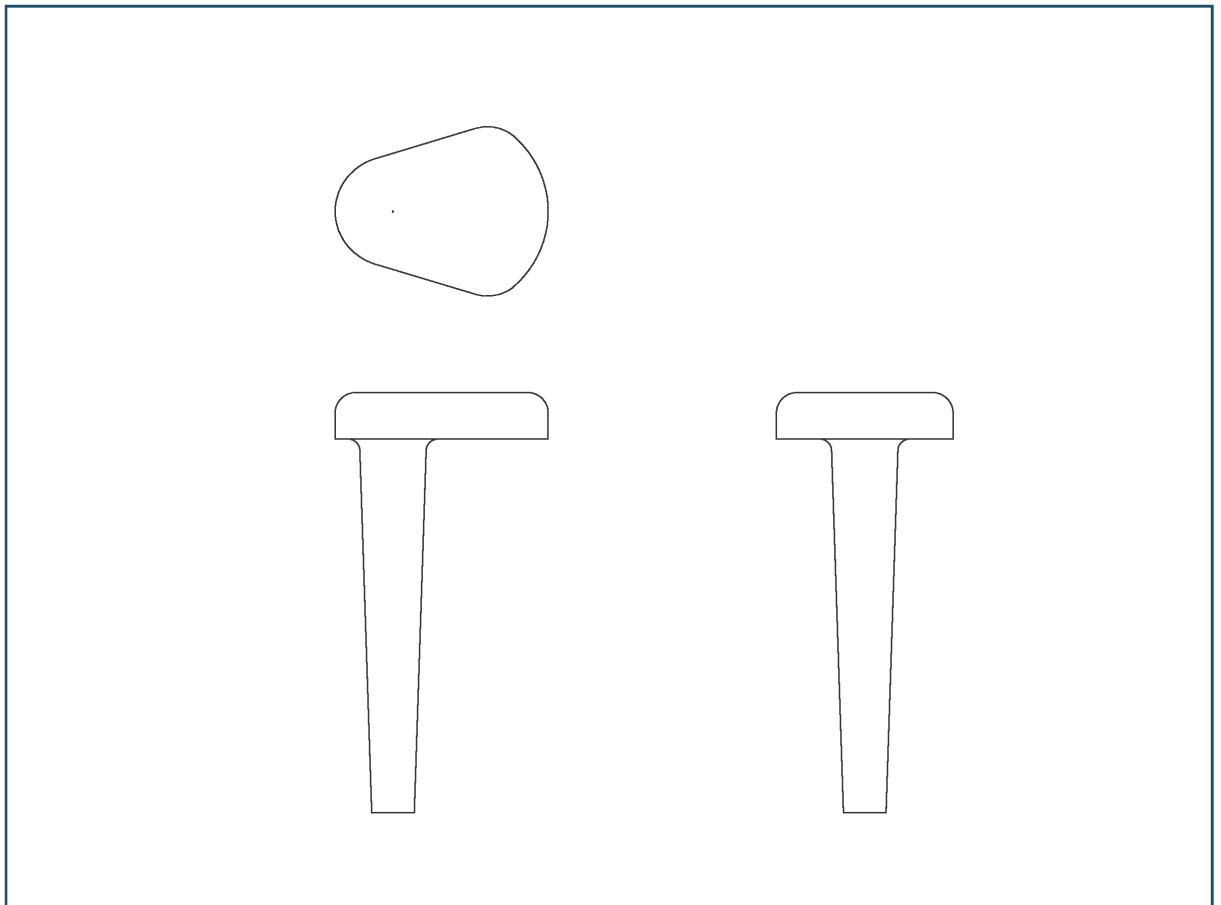


Dimensions

Total length (mm):
5.00 - 5.50 - 6.00
6.50 - 7.00 - 7.50 - 8.00
Flange maximum length: 2.00 mm
Flange thickness: 0.80 mm
Columella distal end diameter: 0.65 mm

Packaging

Dispenser box containing 6 sterile, individually packaged implants.

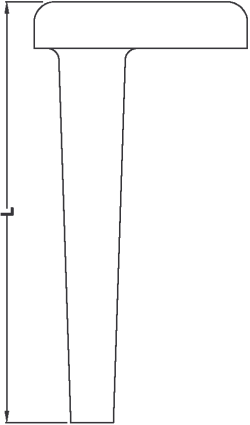	Total length (L) mm	Code
	5.00	TPL 07.05.500
	5.50	TPL 07.05.550
	6.00	TPL 07.05.600
	6.50	TPL 07.05.650
	7.00	TPL 07.05.700
	7.50	TPL 07.05.750
	8.00	TPL 07.05.800

